

THZ-D


THz Detectors for use with our universal monitors



AVAILABLE MODELS



THZ12D-3S-VP
(3W-Volume Absorber)

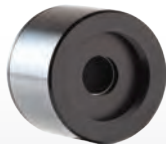

THZ9D-20mS-BL
(9 mm-Organic Black)

KEY FEATURES

1. **RELATIVE MEASUREMENTS FROM 0.1 TO 30 THz**
Broadband, room temperature operation, easier to use and less expensive than a Goly cell
2. **FLAT SPECTRAL RESPONSE**
Get the best precision across the entire wavelength range
3. **MEASURE HIGHER POWERS**
Up to 3 W of continuous power with the THZ12D model, the highest in our terahertz range of products
4. **USE WITH A UNIVERSAL MONITOR**
No need for an exclusive monitor. These unique THz detectors work with our standard universal monitors:
 - MAESTRO
 - M-LINK
5. **LARGE APERTURES**
Models range from 9 to 12 mm \varnothing aperture
6. **CALIBRATED AT 10.6 μ m**
THZ-D detectors are calibrated at a single wavelength (10.6 μ m) and include typical wavelength correction data from 10.6 to 440 μ m. They are used for relative measurements outside that range

ACCESSORIES



Stand with Steel Post
(Model Number: 200160)



Extension Cables
(4, 15, 20 or 25 m)



SDC-500 Digital
Optical Chopper



Pelican Carrying Case

SEE ALSO

HOW IT WORKS	14
TECHNICAL DRAWINGS	126
ABSORPTION CURVES	128
COMPATIBLE MONITORS	
MAESTRO	20
M-LINK	32
LIST OF ALL ACCESSORIES	190



THZ-D

SPECIFICATIONS

	THZ9D- 20mS-BL	THZ12D-3S-VP
MAX AVERAGE POWER	25 mW	3 W
EFFECTIVE APERTURE	9 mm Ø	12 mm Ø
COMPATIBLE MONITORS	MAESTRO, M-LINK & APM	MAESTRO & M-LINK
MEASUREMENT CAPABILITY		
Spectral Range ^a		
Frequency	0.1 - 30 THz	0.1 - 30 THz
Wavelength	3000 – 10 µm	3000 – 10 µm
Maximum Average Power		
with MAESTRO	20 mW	3 W
with M-LINK	25 mW	3 W
Noise Equivalent Power ^b	300 nW	0.5 µW
Minimum Measurable Power ^c	N/A	50 - 100 µW
Thermal Drift ^d	N/A	12 µW/°C
Rise Time (nominal) ^d	<0.2 sec	3 sec
Sensitivity (typ into 100 kΩ load) ^e	120 V/W	200 mV/W
Minimum Repetition Rate ^d	1000 Hz	7 Hz
Chopping Frequency	10 Hz (required)	N/A
Calibration Uncertainty ^f	±5.0 % @ 10.6 µm; ±15 % @ 10.6 - 440 µm ^a	±3.0 % @ 10.6 µm; ±15 % @ 10.6 - 440 µm ^a
Repeatability	±0.5 %	±0.5 %
DAMAGE THRESHOLDS		
Maximum Average Power Density ^g	50 mW/cm ²	30 W/cm ²
Maximum Energy Density	<0.1 J/cm ²	<1 J/cm ²
PHYSICAL CHARACTERISTICS		
Effective Aperture	9 mm Ø	12 mm Ø
Absorber (High Damage Threshold)	BL (Black Absorber)	VP (Volume Absorber)
Dimensions	38.1Ø x 26.2 mm	73H x 73W x 28D mm (80D mm with tube)
Weight (head only)	91 g	320 g
ORDERING INFORMATION		
Product Name	THZ9D-20mS-BL	THZ12D-3S-VP-DO
Product Number		202229

a. From 10 to 440 µm, spectrometer measurement with multiple laser references validation.
From 440 to 600 µm, spectrometer measurement only.
From 600 to 3000 µm, relative measurement only.
This spectral range is subject to change.

b. Nominal value, actual value depends on electrical noise in the measurement system.

c. Actual value depends on ambient conditions and the measurement system.

d. With Gentec-EO monitors.

e. Maximum output voltage = sensitivity x maximum power.

f. Including linearity with power.

g. At 1064 nm, 1 W CW.

Specifications are subject to change without notice